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NOT / NOTES

English

Usage Instructions 20
Meaning of Symbols 20
Responsibility 20
CE Declaration of Conformity 20
Intended Use 20
Unauthorized Use 21
Safety Warnings 21
Product-Specific Safety Warnings 22
Opening the Box 23
Features/Specifications 24
Performance Index 24
To Install or Replace Battery 24
Key Illustration 25
LCD Display 25
To Set Weight Unit 26
Optional Normal Weighing Mode (bypassing parameter setting) 26
Fat/Hydration/Muscle/Bone/Calorie Measuring Mode 27
Parameter Setting 28
Warning Indications 30
About Fat 31
About Calorie 32
Advise for use & care 32
Device Storage 32
Recycling 32
Troubleshooting 33
Shipping 33

Dear Customer,

We hope you enjoy using your product. Thank you for choosing **Fakir Hausgeräte**, a leader in innovative technologies. Please read this user manual carefully before using the product and keep it for future reference.

Usage Instructions

These usage instructions are prepared to ensure that the device is used properly and safely. For your safety and the safety of others, please read the user manual before operating the device.

Meaning of Symbols

In this user manual, you will see the following symbols:

WARNING!

This symbol indicates hazards that could lead to injury or death.

CAUTION!

This symbol indicates hazards that could cause damage to your product or malfunction.

NOTE!

Indicates important and useful information in the user manual.

Responsibility

Fakir Elektrikli Ev Aletleri San. ve Tic. A.Ş. accepts no responsibility for any damage or repairs made outside Fakir Authorized Service Centers or for improper use of the device resulting from non-compliance with the instructions in this user manual.

CE Declaration of Conformity

This device complies fully with the European Union's 2014/30/EU Electromagnetic Compatibility and 2014/35/EU Low Voltage Directives. The product bears the CE mark on the device type label.

Intended Use

The product model is designed for home use. It is not suitable for workplace or industrial use. Any malfunction resulting from such use will not be covered by our authorized service.

Unauthorized Use

In the following cases of unauthorized use, the device may be damaged or cause injuries:

- **Restricted Capabilities and Supervision:** This device should not be used by persons (including children) with limited physical, sensory, or mental capabilities or lack of experience and/or knowledge, unless they are supervised by a person responsible for their safety or given instructions regarding the safe use of the device. Children should be supervised to ensure they do not play with the device.
- **Use with Limited Abilities:** The device can be used by children aged 8 years and above and persons with limited physical, sensory, or mental capabilities or lack of experience and knowledge only if they are supervised, given instructions on safe use, and understand the potential hazards.
- **Packaging Materials:** Keep materials like bags included in the package away from children. There is a risk of suffocation!
- **Cleaning and Maintenance:** Children under 8 years should not clean or maintain the device unless they are supervised.
- **Device and Cable Safety:** Keep the device and its cable out of the reach of children under 8 years old. Children should be supervised to ensure they do not play with the device. Electrical devices are not toys for children. Therefore, use and store the device out of the reach of children. Children may not be aware of the dangers associated with the use of an electrical device. Ensure the power cord does not hang down from the device.

Safety Warnings

When using an electrical device, always follow the precautions below to prevent fire, electric shock, and personal injury:

- **Read the User Manual:** Carefully read the information in the user manual and keep it for future reference. If you give the device to a third party, be sure to include the manual.
- **First Use:** When you remove the product from its packaging, make sure all accessories and attachments are complete. If any parts are missing,

contact your nearest Fakir Authorized Service Center. Keep the protective bags of the accessories and attachments away from children and infants to prevent the risk of suffocation.

- **Check the Condition of the Device:** If you suspect that there is any issue with the device, its casing, cable, or plug, do not operate the device and contact the nearest Fakir Authorized Service Center.
- **Cable and Plug Safety:** Do not use the device if the power cord or plug is faulty or damaged. Do not use the device with an extension cord and ensure that the power cord does not come into contact with heat sources or sharp edges. When unplugging, always pull the plug itself, not the cord. Do not pull on the cord or lift the device by the cord.
- **Child Safety:** Never leave the machine unattended while it is plugged in. Do not allow children to play with the device unsupervised. Store the device's power cord out of the reach of children. Electrical devices are not toys for children.
- **Use of the Device:** Do not use the device outdoors. Your machine may be adversely affected by dirt or rain. Do not immerse the machine in water and protect it from rain or moisture. The device operates with electrical power, which poses an electric shock hazard. Therefore, do not touch the plug with wet hands. Do not use the device to clean flammable or explosive materials (e.g., gasoline or alcohol), burning or smoldering objects (e.g., cigarettes, matches, or hot ashes).
- **Maintenance and Cleaning:** Before cleaning and maintenance, unplug the device and allow it to cool down.
- **Usage Restrictions:** Use the electrical device only as intended. Using the device contrary to the instructions may shorten its lifespan and pose safety risks.

Product-Specific Safety Warnings

User-friendly Tips

1. During Measuring

A. The device is designed to allow auto step on. Please always ignore the first reading and only start taking effective record from the second time weighing.

B. Always remove your shoes and socks , wipe clean your feet before using.

C. Step on the weighing platform gently.
D. It is advisable to take measurement at the same time of the day.
E. Measured result could be misleading after intensive exercise, excessive dieting or under extreme dehydration condition.
F. Always weigh/use the scale on a hard and flat surface. When you are measuring, please keep still.
G. The measured data of the following people may have deviation:
- Children under 10 or adult over 100 (can only use the device under normal weighing mode)
- Adult over 70 Body building athletes or other occupational athletes.

H. People with the following conditions should not use the device:
- Pregnant women
- Those with symptom of edema
- Those on dialysis treatment
- Those who use heart pacemaker or with other implanted medical device.

2. Safety Warning

A. The Fat% Hydration% content evaluated by the Fat Monitor Scale is only for reference (not for medical purpose). If your Fat% Hydration% content is over/below normal level, please consult your doctor for more advice.

B. Keep away from strong electro-magnetic field when using the scale.

C. It could be slippery if you weigh on a wet surface. Moreover, if the weighing platform surface is glazed, never step on the edge of the scale and do keep balance during your weighing.

Carefully read the safety rules and usage instructions in this user manual.

Keep this user manual for future reference.

Opening the Box

- Before first use, remove the device and accessories from the box and check for damage or defects. If the device is damaged or defective, do not use it and contact **Fakir Authorized Service**.
- Repairs should only be performed by Fakir Authorized Service Centers.
- Repairs by unauthorized persons may be harmful to the user.

⚠ WARNING!

Do not use a faulty or defective product.

Features/Specifications

1. Using high precision strain gauge sensor.
2. Multifunction: test body fat %, body hydration %,body muscle%, body bone% and suggest calorie intake requirement.
3. Built-in memory for 12 different users.
4. Optional normal weighing mode (bypassing parameter setting): Auto step on weighing and consecutive add-on weightings.
5. Low battery and overload indication.

NOTE!


The values indicated on the device or in other printed documents provided with the product are obtained under laboratory conditions according to relevant standards. These values may change depending on the product's use and environmental conditions.

Performance Index

Index / Item	Weighing Capacity	Fat%	Hydration%	Muscle%	Bone%	Calorie
Index Range		4.0% ~60.0%	27.5% ~66.0%	20.0% ~56.0%	2.0% ~20.0%	
Division	0.1kg 0.2lb	0.1%	0.1%	0.1%	0.1%	1KCAL

To Install or Replace Battery

- Remove isolation sheet before use.
- Open the battery case cover at the bottom of scale.
- Press battery gently, the battery will pop up automatically.
- Install new battery by connecting one side of battery to the positive plate first, and then pressing down the other side. Please observe the polarity.

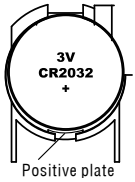
Key Illustration

SET (ON/SET) : Turn on the scale and prompt to next parameter. (Please refer to the Setting Personal Parameter overleaf)

▲: Advance the value by one interval each press. Continuous pressing will speed up the advance of value.

▼:Decrease the value by one interval each press. Continuous pressing will speed up the decrease of value.

LCD Display

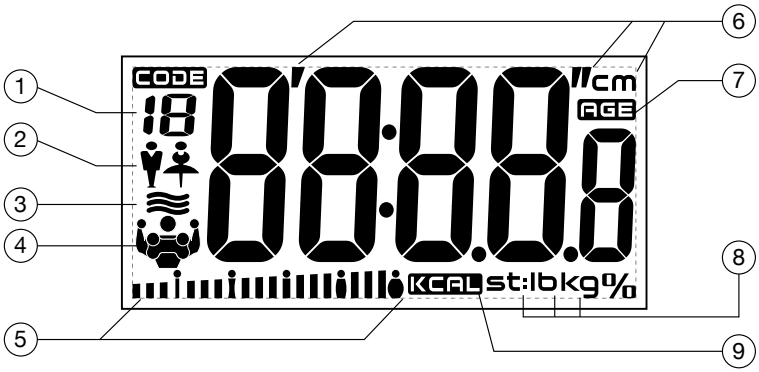
1. User memory code
2. Gender
3. Hydration%
4. Fat%/Muscle%Bone%
5. Fat status indication
6. Height Unit
7. Age
8. Weight Unit
9. Calorie Unit

Here is the measurement assets and the sorting the names of the measurement percentages.


While the measurement;

- First measurement is the weight (Kg).
- Second measurement is %Fat ratio and status of the fatness in below. Third measurement is % Water ratio. Fourth measurement is % Muscle ratio, fifth measurement is % Bone ratio and the sixth measurement is daily recommended calorie (kcal).

To Set Weight Unit

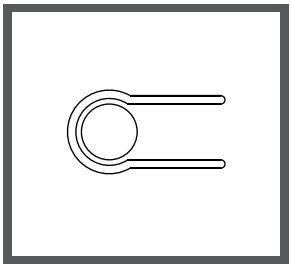KG-Kilogram
LB-Pound
ST-Stone

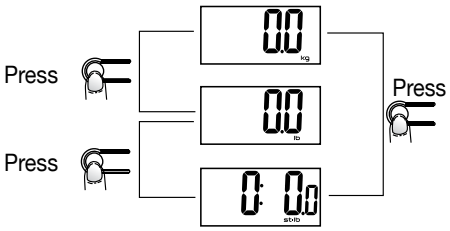At the bottom of scale,you can find a unit conversion button as shown in picture.


- If your scale is available (or applicable) in kg/lb/st, you can set to the one you prefer by pressing the unit conversion button. Press the unit conversion button when scale is power off, LCD shows current weight unit. Press the button again to convert weight unit.
- After selecting the prefered unit, the scale will be auto-off after 3 seconds. Selected weight unit will be activated when the scale is turned on again.
- When setting lb/st as weight unit, the default height unit is feet/inch.
- When setting kg as weight unit, the default height unit is cm.

NOTE!

The unit system of scale depends on sales destination or the requirement of distributor.

Optional Normal Weighing Mode (bypassing parameter setting)

1. Auto Step On Weighing

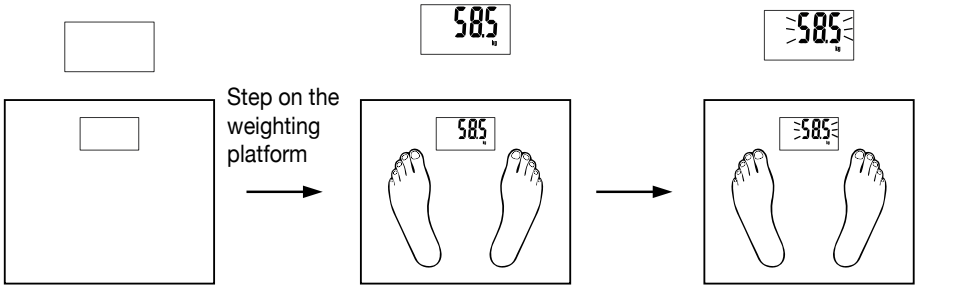Install battery (please refer to Page 2), and place the scale on flat floor.


Display your weight.

Digits flash, stabilized and locked.

2. Start Measurement

After your weight value is stably displayed on the LCD screen, you can perform measurements with additional weights. However, for a new measurement, the added weight must be more than 2 kg.

Example:If you want to obtain your baby's weight, you can first weigh yourself, and then, without stepping off the scale, hold your baby in your arms. Then you will get the total weight reading. By subtracting your previous reading from the new one, you obtain your baby's weight as the difference.

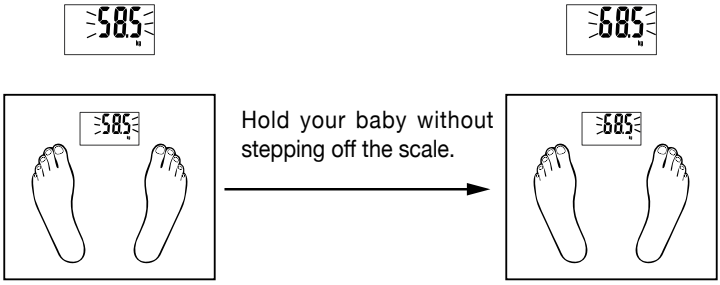Digits flash, stabilized and locked. Display total weight.

Fat/Hydration/Muscle/Bone/Calorie Measuring Mode

1. Setting personal parameters

At the first time when you install battery and switch on, all default parameter are set as follow:


Memory	Gender	Height	Age
1	Male	165cm/5'05.0"	25

Parameters can be set in the following range:

Memory	Gender	Height	Age
1-12	Male/Female	100-250cm 3'03.5"-8'02.5"	10-100

Parameter Setting

Example: The 3rd user, female, 168cm, 28 years old, can set her parameters as follow:

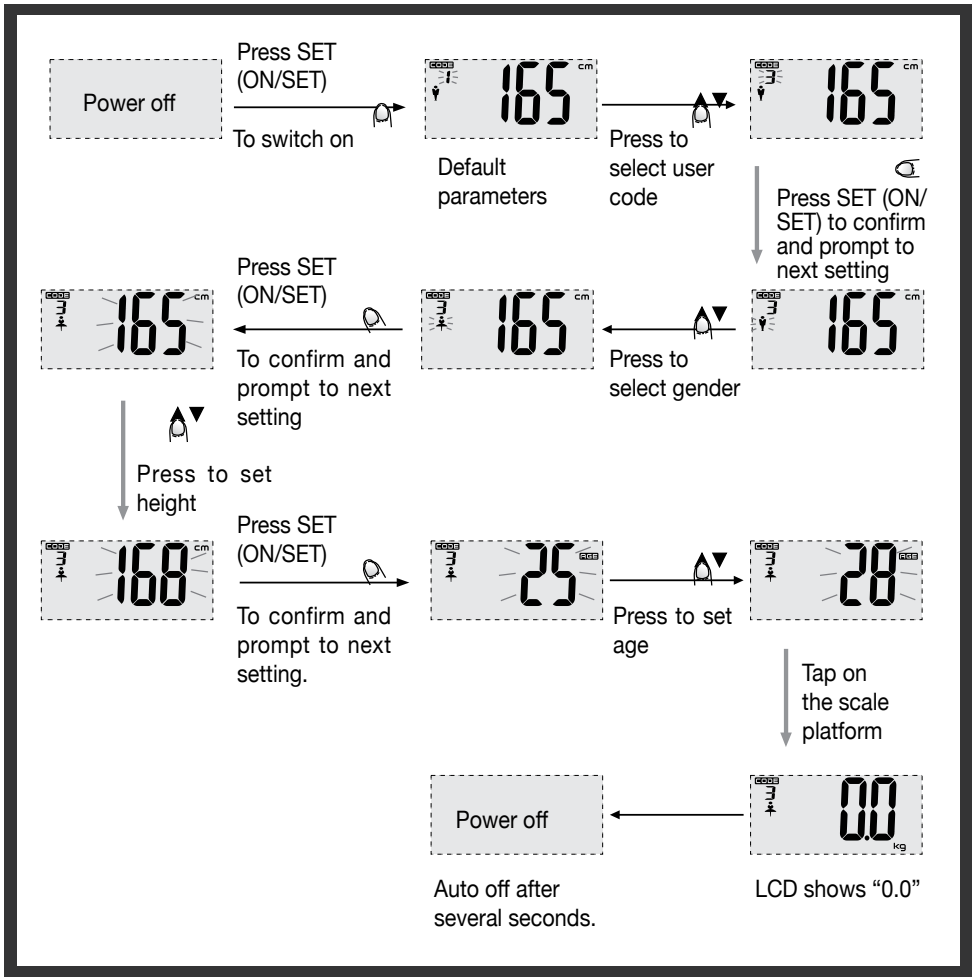NOTE!

During parameter setting, if idle for 6 seconds, the parameters already set will be memorized automatically and LCD shows “0”, meaning body fat measuring will be started.


NOTE!

If LCD prompt to show “0” during parameter setting, press SET(ON/SET) for 3 seconds to return setting mode, then you can set user code,gender, height and age.


2. To Start Measuring

- Make sure your feet are touching with metal electrodes on the scale platform, or body fat can not be measured.
- To ensure accuracy, please remove your shoes & socks before stepping on.

Example: The 3rd user, female, 168cm, 28 years old can measure as follow.

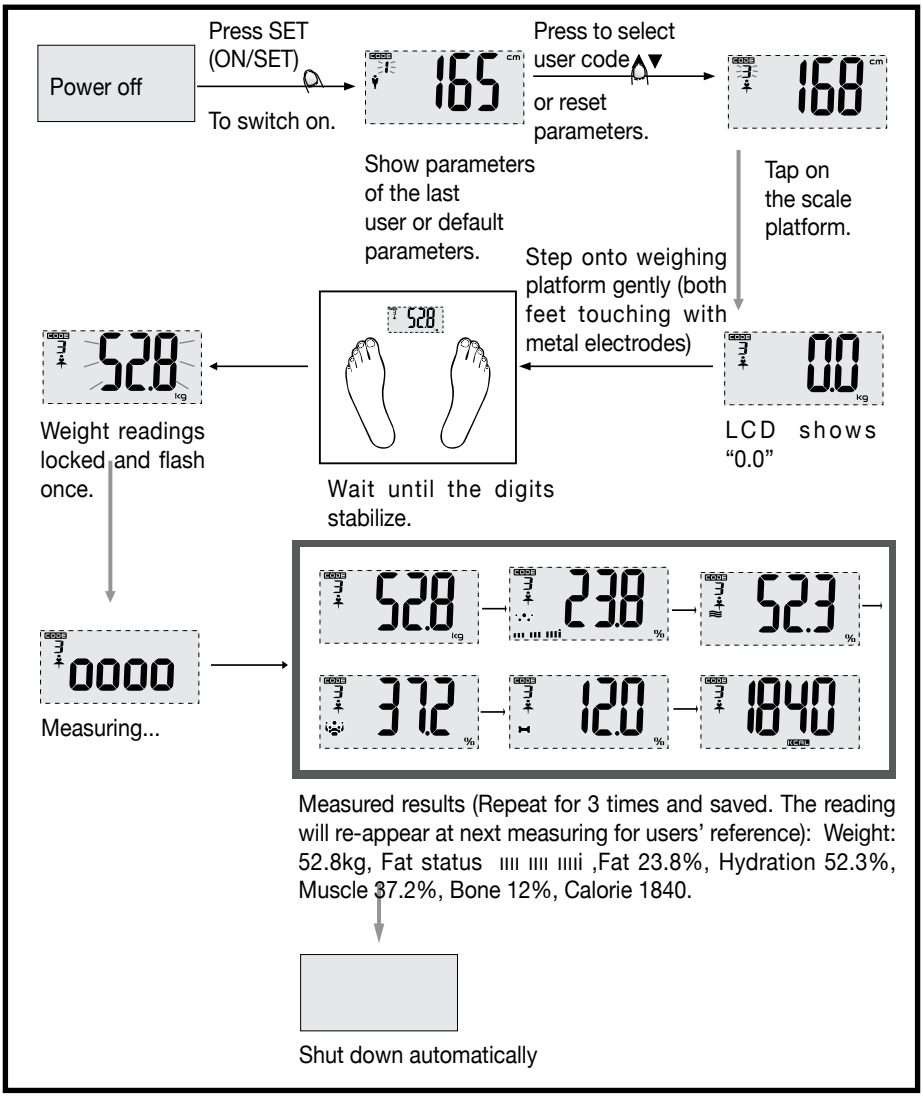
Warning Indications

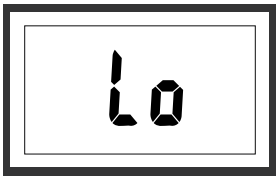Low Battery Indication

The battery power is running low, please replace with a new battery.

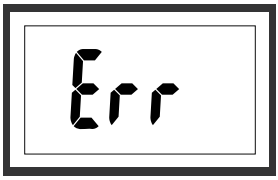
Over-load Indication

The weighing subject on the platform exceeds the maximum capability of scale. Please step off to avoid damage.


Low Fat% Indication

The Fat% is too low. Please lead a more nutritional diet and take good care of yourself.

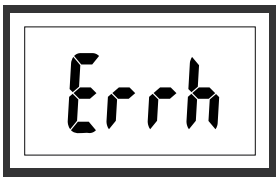High Fat% Indication

The Fat% is too high. Please watch your diet and do more exercise.

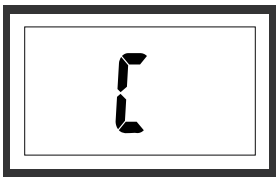Re-test

Error exists, please re-test to get the correct result.

About Fat

1. About Fat%

- Fat% is an index % of fat content in human body.
- The fat content in human body is of a certain level, overabundance of fat will bring bad effects on metabolism, thereby, various illnesses may invade. Here goes an old saying, “Obesity is hotbed of illness”.

2. Working Principles

As electrical signal of certain frequency transmits through human body, the impedance of fat is stronger than that of muscle or other human tissues. Through the use of biological resistance analysis, this product will pass a safe frequency of low voltage electrical current through human body, thus, the Fat% will be reflected according to the resistance differences when the electrical signal transmitting through human body.

3. Fat/Hydration Content Fitness Assessment Chart

After your weighing, you can use the following Assessment Chart for reference:

Age	Fat(F)	Hydration(F)	Fat(M)	Hydration(M)	Status	Key Line Indication
≤30	4.0-16.0%	66.0-57.8%	4.0-11.0%	66.0-61.2%	Too Lean	ıiiii
	16.1-20.5%	57.7-54.7%	11.1-15.5%	61.1-58.1%	Lean	iii iiii
	20.6-25.0%	54.6-51.6%	15.6-20.0%	58.0-55.0%	Normal	iii iii iiii
	25.1-30.5%	51.5-47.8%	20.1-24.5%	54.9-51.9%	Fat	iii iii iii iiii
	30.6-60.0%	47.7-27.5%	24.6-60.0%	51.8-27.5%	Too Fat	iii iii iii iii iiii
≥30	4.0-20.0%	66.0-55.0%	4.0-15.0%	66.0-58.4%	Too Lean	ıiiii
	20.1-25.0%	54.9-51.6%	15.1-19.5%	58.3-55.3%	Lean	iii iiii
	25.1-30.0%	51.5-48.1%	19.6-24.0%	55.2-52.3%	Normal	iii iii iiii
	30.1-35.0%	48.0-44.7%	24.1-28.5%	52.2-49.2%	Fat	iii iii iii iiii
	35.1-60.0%	44.6-27.5%	28.6-60.0%	49.1-27.5%	Too Fat	iii iii iii iii iiii

The above information is only for reference

About Calorie

1. **Calorie & Energy Consumption**
 Calorie: Energy Unit
 The data reading indicates the energy you need to sustain the daily metabolism and ordinary activities required for your weight, height, age and sex.
2. **Calorie & Losing Weight**
 Control the calorie intake combined with proper exercise is an effective way to lose weight. This is a very simple theory. As the daily energy intake can not support the energy consumption of the human body, the saccharide and fat stored inside will decompose and provide energy to human body. That leads to weight loss.

Advise for use & care

- The platform will be slippery when it is wet. Keep it dry!
- Keep standing still during the course of weighing.
- DO NOT strike, shock or drop the scale.
- Clean the scale with a damp cloth and avoid water getting inside the scale. DO NOT use chemical/ abrasive cleaning agents.
- Keep the scale in cool place with dry air.
- Always keep your scale in horizontal position.
- If failed to turn on the scale, check whether battery is installed, or battery power is low. Insert or replace new battery.
- If there is error on display or can not turn off the scale for a long time, please take off the battery for about 3 seconds, then install it again to re move soft-fault. If you can not resolve the problem, please contact your dealer for consultation.

Device Storage

- If you do not plan to use the appliance for a long while, store it in a cool and dry place.
- Ensure that the appliance is unplugged and its parts are completely dry.
- Ensure that the appliance is kept out of reach of children.

Recycling

When your product reaches the end of its life, follow these steps for proper recycling:

- Cut the device's cable to make it unusable.
- Dispose of the device in designated waste bins for such devices according to the applicable laws in your country.

- Electrical waste should not be disposed of with household waste.
- Throw the device into bins designated for electrical waste.
- The device's packaging is made from recyclable materials.
- Dispose of packaging materials in appropriate recycling bins.
- The device complies with the 2
- The device complies with the WEEE (Waste Electrical and Electronic Equipment) Directive.

NOTE!

- Check if you can fix the problem or defect yourself before calling customer service. For this, please refer to the "**Troubleshooting**" section. If you cannot resolve the issue or defect, call customer service.
- Before calling, please have the **model name** and product **serial number** ready. You can find this information on the **product's type label**.

 Help Hotline: (+90 212) 444 0 241

Troubleshooting

Your appliance has undergone quality control to prevent any malfunctions. However, if it does not operate for any reason, please check the warning indicators before contacting an authorized **Fakir Service Center.**

Shipping

To prevent damage to the device, follow these steps to transport your device safely:

- Transport the device in its original packaging.
- If the original packaging is not available, use a suitable and soft packaging that protects the device.

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Scan the QR code for the user manual in different languages.

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